



Thermal Innovations and Products

# Single TEC Control

► [Datasheet](#)

# Single TEC Control

## Product Description

The INHECO Single TEC Control (STC) is a compact plug & play temperature and RPM controller for the use with INHECO's CPAC, HeatPAC and temperature controlled shaking devices. The controller box includes a temperature controller with high accuracy, a device connector slot, a touch display, a wide-range power supply for worldwide use, an USB interface, a rugged housing and a RPM controller for shaking operations.

The STC can be individually equipped with a controller slot for either a CPAC, HeatPAC or temperature controlled shaker device. The device slot can be exchanged easily, giving maximum flexibility. The STC can be either operated manually using the touch screen, or through an USB interface. The soft- and firmware of the controller allows programming of complex temperature and shaking sequences.

## Application Areas

The STC is mostly used together with the principal INHECO heating/cooling/shaking units in/on robotic platforms and systems in Lab Automation. The purpose is to form temperature controlled areas with direct and easy access by robotic grippers and pipetting heads. The STC is designed for the use in fundamental research, drug discovery, clinical research, diagnostics and other Life Science applications. The STC is also prepared for use in customer IVD applications.

The STC is used with the following INHECO heating/cooling units:

- CPAC Microplate
- CPAC Ultraflat
- HeatPAC
- Variomag Teleshake
- Variomag Thermoshake

For further information on these units, see the relevant data sheets, which can be downloaded from our website.

## Variants

The INHECO STC is specially designed for controlling a single temperature controlled position. For multiple controlled positions we refer to the datasheet Multi TEC Control (MTC), which can be downloaded from our website.



8900031: INHECO STC



8900031: Backside

Product	Article No.	Dimension W x D x H (mm)	Weight (kg)
Single TEC Control STC	8900031	147 x 224 x 177	3,3 kg

## Verification Plate

INHECO offers a compact and precise measurement plate (IMP) in SBS format for verification of temperature and/or shaking performance. The IMP can be used on all INHECO LabAutomation heating, cooling, shaking and incubating devices. We refer to the IMP datasheet for more information.



7901000: INHECO Measurement Plate


## Main Specifications

Specification	8900031
Input	100-240Vac, 50/60 Hz
Output	1x 24Vdc, 4,5A
Accuracy	±0,5°C
Output CPAC Slot	12-14V, 4,5A
Output Shaker Slot	24V 4,5A
Ambient temperature at max. power	15°C to 32°C
Interface	USB (1.1 or 2.0)



8900031: Easy access to device slot

## Drawings

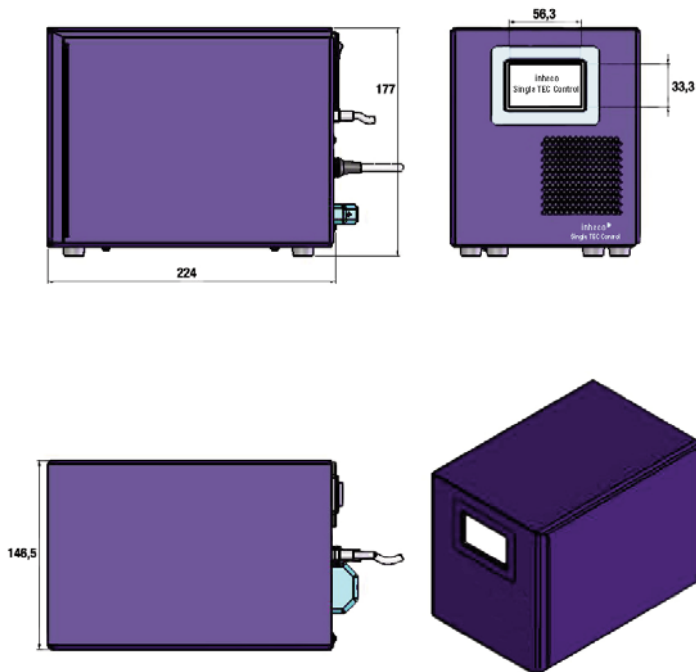

## Accessories & Related Products

Product	Article No.
Slot CPAC / CPAC HT	2400128
Slot Shaker / CPAC HT 2-TEC / HeatPAC	2400125



2400125 Device Slot Shaker

**INHECO**  
Industrial Heating & Cooling GmbH

Fraunhoferstrasse 11  
82152 Martinsried  
Germany

Tel.: +49 (0)89 899593-101  
Fax: +49 (0)89 899593-149

E-mail: [sales@inheco.com](mailto:sales@inheco.com)  
[www.inheco.com](http://www.inheco.com)

This Datasheet can also be downloaded from our website. For further information, please contact our sales department.

